

GETINGE
GETINGE GROUP

**GETINGE WD14 TABLO
A FAST, HIGH-THROUGHPUT TABLETOP
WASHER-DISINFECTOR**



Always with you

IT'S ALL ABOUT TIME

Small clinics and dental practices have big demands – on being able to reuse instruments and circulate items quickly. And of course with complete reliability, safety and cost-efficiency.

We have the perfect solution: Getinge WD14 Tablo – the new generation of tabletop washer-disinfectors.



Getinge has over 100 years of knowledge in infection control and contamination prevention. Our experience comes from installations in everything from sterilization centers and huge hospitals to small health-care units and private medical/dental facilities. We apply the experience gained here in the development of each new product.

Our focus is on the needs of our customers and how we can help to improve care, cut costs, reduce strain, ensure hygiene and enhance people's quality of life.





90° POWER CLEAN 20:00
Cycle Clean Start



GETINGE
Washing Machine



Clean 20:00 20:00
90° POWER CLEAN



GETINGE
Washing Machine



SUPREME PERFORMANCE IN A SMALL COSTUME

A tabletop washer-disinfector with modern design, multifunctional washing features and compliant with all relevant norms and standards.

COMPACT

The compact size, 500x527x657 mm (WxHxD), makes Tablo easy to place, operate and maintain. It has also a unique chamber design with large loading capacity and long spray wing (450 mm) for perfect results.

Placing the washer-disinfector on a standard worktop gives you a more ergonomical working height.

FAST

The speed of the process can play a key role in your overall economy, as it determines when instruments are available again – and how many instruments you need to have in stock. Tablo can cut your process time to just 40 minutes with hot water connection and 3-phase electrical current.

An integrated ventilation system removes the residual moisture to dry the goods. It is possible to interrupt the ventilation phase after 5 minutes if needed. It also cuts your water and energy consumption, adding to cost-efficiency. Even the installation is fast.

RELIABLE

Tablo is a fully automatic washer-disinfector with independent temperature monitoring that supervises the complete process and, in case of malfunction, gives an alarm to the user.

For added safety, an electrical locking system prevents opening the door during a process. And there's a control system for trouble-shooting.

You only have to load the machine, set the program and press the button.

WELL-TESTED

Getinge's range of detergents has been developed, tested and validated on Getinge washers to consistently meet EN ISO 15883 standards for cleaning efficiency of goods commonly found in CSSD's, clinics and practices to protect users and patients.



OPTIONS



DOSING PUMPS

Tablo is equipped with 2 dosing pumps with flow control for detergent as an option. Two extra dosing pumps are also available for combining with enzyme, lubrication and neutralizer.



PRINTER

The printer provides a detailed report with all information about the completed cycle.



NETCOM

Tablo has an RS232 port to connect to the USB storage solution for process documentation. The USB storage solution is available as an option.

A VARIETY OF ACCESSORIES

Getinge offers a complete range of accessories specially designed to assure efficient cleaning.

DENTAL



BASIC CART w. BLIND DOCKING

Blind docking comes as standard. When choosing MIS or Luer Bar the cart will be delivered without blind docking.



STAND FOR INSTRUMENTS

The stand for instruments can hold up to 16 dental instruments.
Art. no: 4000004-01




DENTAL KIT

With 1 tray rack for 9 dental trays, 1 SPRI III instrument tray and 1 fine mesh basket with lid for small parts. (Trays in pictures not included)
Art. no: 502795400



DENTAL BAR

Equipped with 10 injectors for hand pieces. 5 for hand pieces with an approx. outer diameter of 16 mm and 5 for hand pieces with an approx. outer diameter of 20 mm. Fine filter is included. (Hand pieces in pictures not included)
Art. no: 501261200

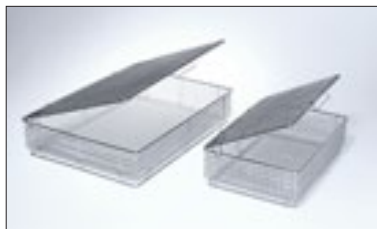

INSTRUMENT TRAY

Art. no: 466 30 91-01 (250x170x70 mm)
Art. no: 466 30 92-01 (340x250x70 mm)
Lid:
Art. no: 498 63 74-00 (250x170x70 mm)
Art. no: 477 34 53-00 (340x250x70 mm)
Lid with handle:
Art. no: 498 81 69-00 (250x170x70 mm)
Art. no: 498 81 69-01 (340x250x70 mm)



FINE-MESH INSTRUMENT TRAY LID

Art. No. 93 28 95 93 28 96 (210x150x40 mm)
Art. No. 93 28 94 476 94 40-00 (310x210x40 mm)
Art. No. 70 53 00 70 53 20 (420x310x40 mm)



AVAILABLE IN TWO COLORS

Getinge WD14 Tablo is available in two colors, white coated polymeric and stainless steel or blue coated polymeric and stainless steel.



Getinge WD14 Tablo

External dimensions	500x527x657 mm (WxHxD)
Internal dimensions	435x300x460 mm (WxHxD)
Weight	35 kg
Volume	63 L
Chamber and inner door materials	AISI 316L
Outer material	Polymeric material (ABS)
Panel material	AISI 304L
Communication ports	RS232 and RS485
Electrical	Single- or 3-phase 50/60 Hz 200V, 208V, 230V or 400V
Heating element	3.6 kW or 1.8 kW
Water supply	Cold or cold/hot
Dosage pumps	2
Circulation pump	115 L/min
Noise	< 60 dB (A)
Standards	EN ISO 15883, MDD, EN IEC 61010-1, EN 61000-3-2
Optional	2 extra dosing pumps with flow control External printer NetCom



COMPLETE SOLUTIONS FOR INFECTION CONTROL

Getinge is the world's leading provider of solutions for effective cleaning, disinfection and sterilization in the healthcare and life science sectors. We are dedicated to helping our customers provide better care at a lower cost. We do this by offering well-thought-through and customized solutions. This means that we are with our customers all the way from architectural planning and education to traceability and support – with complete solutions, long-term commitment and global presence. Getinge – Always with you.

Getinge Infection Control AB

P O Box 69
SE-305 05 Getinge
SWEDEN
Phone: +46 10 335 00 00
info@getinge.com

GETINGE

GETINGE GROUP

www.getinge.com

GETINGE GROUP is a leading global provider of products and systems that contribute to quality enhancement and cost efficiency within healthcare and life sciences. We operate under the three brands of ArjoHuntleigh, GETINGE and MAQUET. **ArjoHuntleigh** focuses on patient mobility and wound management solutions. **GETINGE** provides solutions for infection control within healthcare and contamination prevention within life sciences. **MAQUET** specializes in solutions, therapies and products for surgical interventions, interventional cardiology and intensive care.